

DIN EN ISO 24234 Dentistry –Mercury and Alloys for Dental Amalgam / Determination of Creep , Dimensional Change, Compressive Strength

Description

6.6 Preparation of test specimens to determine compliance with the requirements for creep, dimensional change and compressive strength

6.6.1 Temperature

Prepare specimens at $(23 \pm 2) \text{ }^\circ\text{C}$, unless specified otherwise.

6.6.2 Mixing

Mix a mass of the amalgam sufficient to make a test cylinder (8 ± 1) mm high after condensation into the die shown in Figure 1. If necessary, adjust the height of the cylinder by condensing only a portion of the amalgam mix into the die.

For alloy supplied either as tablets or as a powder in bulk, the alloy mass and the mercury mass to be mixed, the capsule, the pestle (if needed), and any other mixing accessory required shall conform to those

recommended by the manufacturer in accordance with 7.3.2 and 7.2.1 h). Use dental mercury that complies with requirements 4.5 and 4.6 of this International Standard.

For products supplied in predosed capsule form, use as many capsules as needed. If more than one capsule has to be used, mix them simultaneously using equipment of the same type.

Use an amalgamator that complies with ISO 7488. The manufacturer of the amalgam product, in accordance with 7.3.1 b), shall specify the amalgamator's setting and the mixing time.

6.6.3 Apparatus for the preparation of test specimens for determining creep, dimensional change and compressive strength

6.6.3.1 General

Use the apparatus as shown in Figures 1 to 5.


6.6.3.2 Materials and tolerances for construction of the apparatus

Make the holder, the spacers and the cap of cold-rolled or stainless steel. Make the die and the plungers of hardened tool steel or hardened stainless steel. Hone the working surfaces of the die and the plungers to a core roughness depth (Rk) not greater than $6,3 \text{ } \mu\text{m}$ when tested in accordance with ISO 13565-2. Limits of clearance between the die and the plungers shall be F7h7 in accordance with ISO 286-2.

6.6.3.3 Assembly of the apparatus

For production of creep and compressive strength specimens, assemble the holder, spacers Nos. 1 and 2, the die and plunger No. 2 as shown in Figure 1.

Particular measuring instruments used in the dimensional-change test (e.g. interferometers) may require an impression on the end surface of the specimen that is produced by the cap that is shown in Figure 5. For production of dimensional-change specimens, include the cap in the assembly if this is appropriate for the measuring instrument that is to be used. In which case, position the cap on top of the plunger No. 2.

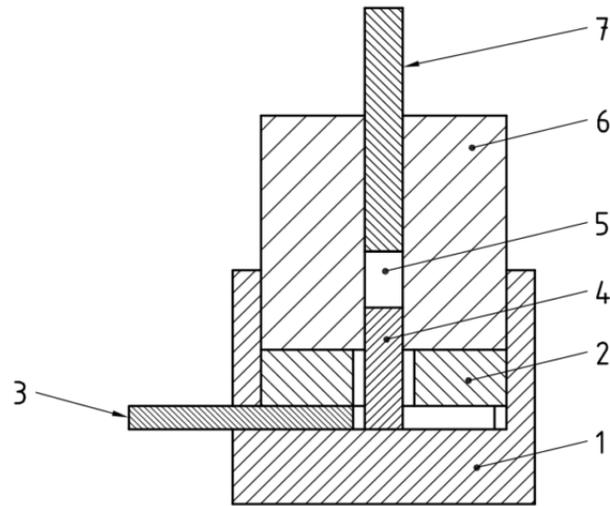


Key

- 1 holder
- 2 spacer No. 1
- 3 spacer No. 2
- 4 plunger No. 2
- 5 test specimen
- 6 die
- 7 plunger No. 1

The dimensions for each of the components are given in the figures that follow.

Figure 1 — Vertical section through the apparatus for making amalgam test specimens, showing the assembled apparatus with a specimen in place

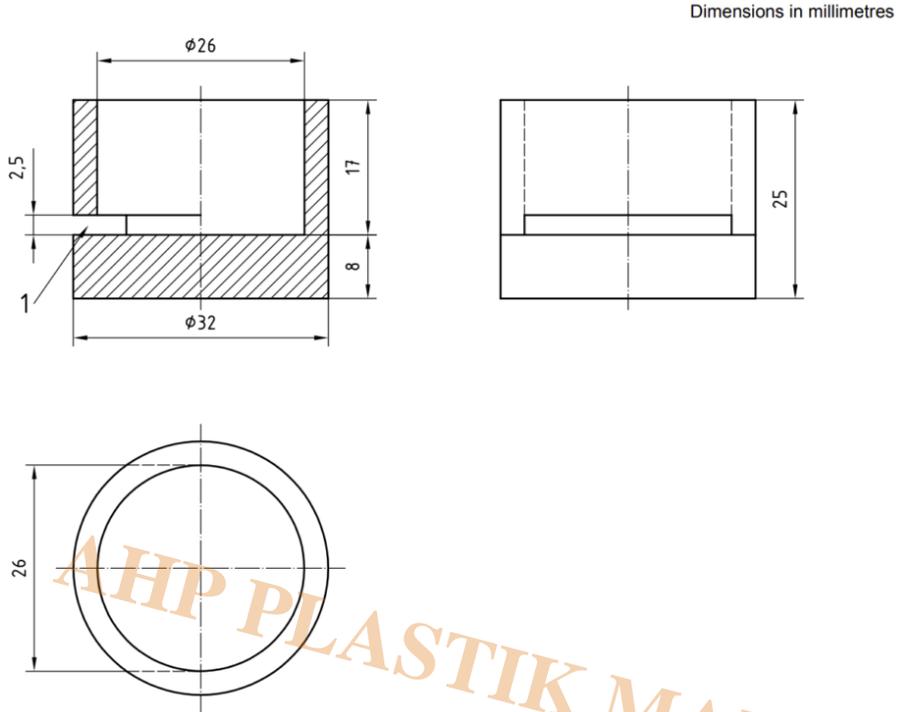

Key

1 slot

Figure 2 — The holder



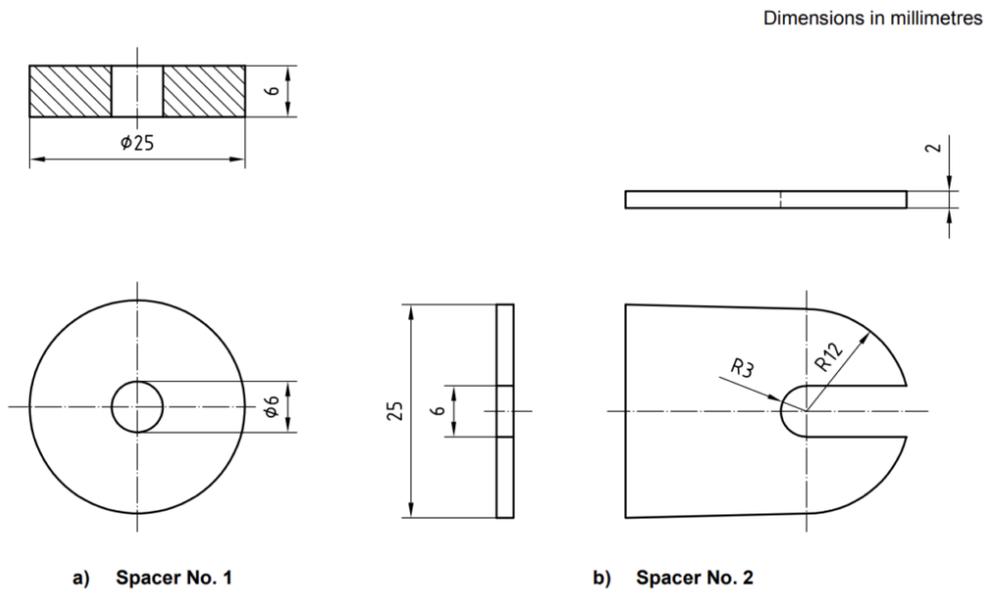


Figure 3 — Spacer No. 1 and Spacer No. 2

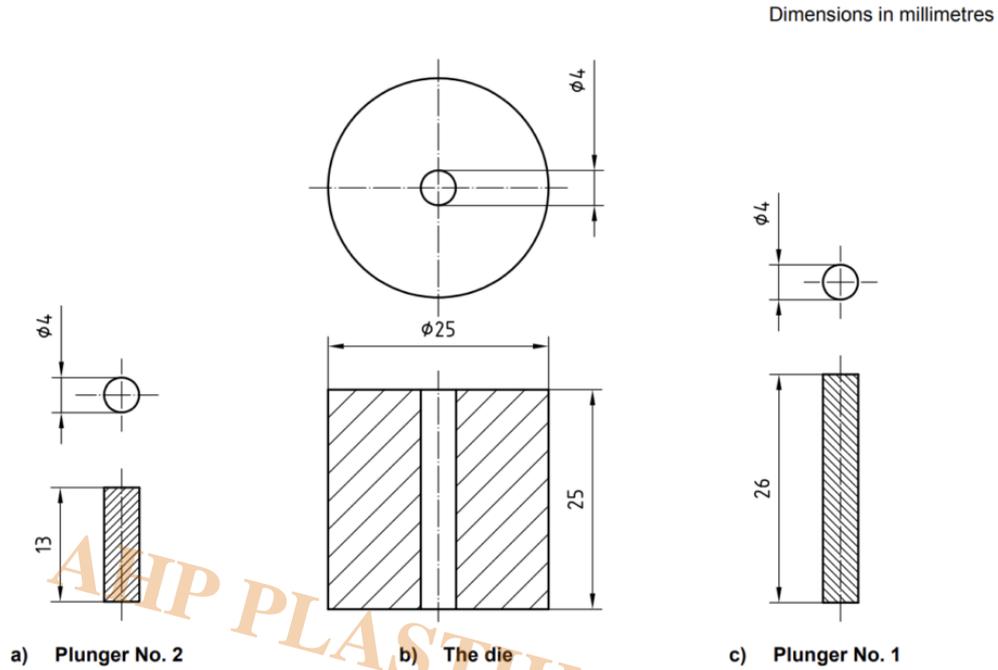

To assist the operator to judge whether the correct quantity of amalgam has been inserted, to give a specimen that is within the permitted range for length, i.e. (8 ± 1) mm, three circumferential datum lines should be scribed at 9 mm, 11 mm and 13 mm from one end of plunger No. 1. This end shall be in contact with the amalgam. Though such datum lines are not mandatory, their use is recommended.

Figure 4 — Plunger No. 2, The die and Plunger No. 1

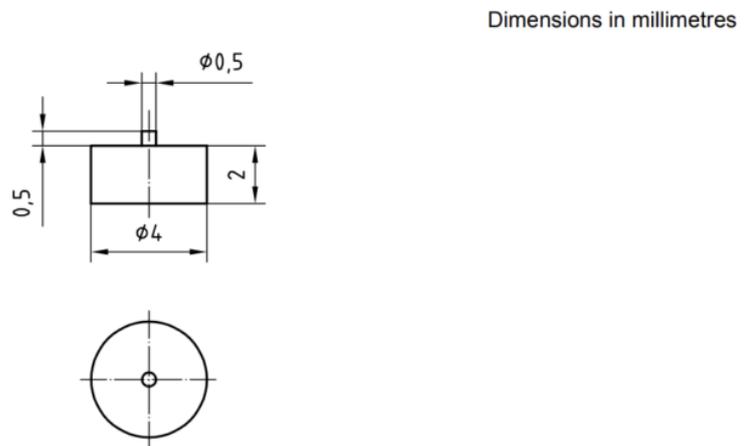

This cap shall be used in conjunction with the apparatus shown in Figure 1 if the measuring instrument requires a small indentation at one end of the specimen.

Figure 5 — Cap used for the production of dimensional-change specimens

Table 3 — Schedule for the preparation of test specimens

Procedure	Time
End of mixing	0
Place mixed mass into the die, insert plunger No. 1 and apply a load to produce a pressure of (14 ± 1) MPa	30
Release load and remove spacer No. 2 at	45
Replace load at	50
Release load at	90
Carefully remove excess mercury and eject specimen at	120

Time in seconds

6.6.4 Condensation

Empty the coherent mass of the mixed amalgam on top of the die cavity and insert immediately, using several thrusts of an amalgam condenser slightly less than 4 mm in diameter. Do not express mercury during insertion.

Insert plunger No. 1 and follow the schedule in Table 3. The specimen shall not be trimmed. After ejection,

transfer the specimen to air maintained at $(37 \pm 1) \text{ }^\circ\text{C}$.

NOTE If plunger No. 1 has circumferential datum lines scribed on its cylindrical surface, they will provide reference points to determine whether the amount of amalgam inserted will give a specimen that is within the permitted range for length. The lines are at 9 mm, 11 mm and 13 mm from the end of the plunger that is in contact with the amalgam.

If the cap (Figure 5) is not present in the assembled apparatus and the 13 mm datum line alone can be seen, the specimen will have a length that is (8 ± 1) mm.

If the cap (Figure 5) is present in the assembled apparatus and both 11 mm and 13 mm datum lines can be seen but the 9 mm line cannot, the specimen will have a length that is (8 ± 1) mm.

6.7 Determination of creep

6.7.1 Preparation of test specimens

Prepare the specimens in accordance with 6.6. Make five specimens. Store these at $(37 \pm 1) \text{ }^\circ\text{C}$ for $(7,0 \pm 0,2)$ days. Prior to testing, grind the surface of both ends of each specimen plane and perpendicular to the axis using wet, coated abrasives that comply with microgrit size P1200, according to ISO 6344-1:1998.

Measure the length of each specimen and record it, to the nearest 0,01 mm, as the original length.

6.7.2 Procedure

Apply a stress of $(36,0 \pm 0,2)$ MPa normally and uniformly over the cylinder ends, continuously for 4 h at a temperature of $(37,0 \pm 0,5) \text{ }^\circ\text{C}$. Record the change in specimen length between $(1,00 \pm 0,05)$ h and

$(4,0 \pm 0,1)$ h to an accuracy of 0,01 mm.

Calculate the creep strain, to the nearest 0,1 %, as follows:

$$\frac{\Delta l}{l_0} \times 100 = \varepsilon_c \quad (3)$$

where

Δl is the change in length occurring between 1 h and 4 h;

l_0 is the original length;

ε_c is the creep strain, as a percentage.

Test three specimens. If all three results meet the requirement in Table 2, it is not necessary to test the other two specimens.

If one of the three specimens fails to meet the requirement in Table 2, test two more specimens.

Test no more than five specimens.

6.8 Determination of dimensional change during hardening

6.8.1 Preparation of test specimens

Prepare the specimens in accordance with 6.6.

6.8.2 Procedure

Place the specimen in the measuring instrument. Do not subject the specimen to a restraint greater than

0,02 N during the test. Maintain the specimen at a temperature of $(37 \pm 1) \text{ }^\circ\text{C}$ during the test.

Measure the dimensional change that occurs between $(5,0 \pm 0,1)$ min and $(24,0 \pm 0,1)$ h from the end of mixing, to an accuracy of $\pm 0,5 \text{ } \mu\text{m}$. At $(24,0 \pm 0,1)$ h measure the specimen length to an accuracy of $\pm 0,01$ mm.

Calculate the dimensional change, to the nearest 0,01 %, as follows:

$$\frac{\Delta l_d}{l_d} \times 100 = \varepsilon_d(\%) \quad (4)$$

where

Δl_d is the dimensional change occurring between 5 min and 24 h;

l_d is the length at 24 h;

ε_d is the dimensional change as a percentage.

6.9 Determination of compressive strength

6.9.1 Preparation of test specimens

Prepare the specimens in accordance with 6.6.

6.9.2 Procedure

Determine the compressive strength by means of a suitable mechanical testing machine. During the test,

maintain the specimen at a temperature of $(23 \pm 2) \text{ }^\circ\text{C}$. Apply the load normally and uniformly over

the cylinder ends at a crosshead speed of $(0,5 \pm 0,1)$ mm/min. For each specimen, record the compressive strength to the nearest 5 MPa.

6.9.3 Compressive strength at 1 h

Determine the compressive strength of five specimens at (60 ± 2) min after mixing.

If only three specimens meet the requirement for compressive strength at 1 h in Table 2, determine the compressive strength of five more specimens. Test no more than ten specimens.

6.9.4 Compressive strength at 24 h

Determine the compressive strength of five specimens at (24 ± 1) h after mixing.

If only three specimens meet the requirement for the compressive strength at 24 h in Table 2, determine the compressive strength of five more specimens. Test no more than ten specimens.



Universal Tensile-Compression Tester According to ISO 24234

- According to ISO 527 , ISO 24234
- Servo controlled
- Ball screw
- Double column
- USB port for computer connection
- Windows based software
- Load resolution 1/10000
- Speed range up to 200 mm/min
- Compression grips for testing of strength of amalgam
- Easy operation and clear visualization (test curves, calculations)
- Easy to change the grips via male-female connection
- Accuracy $\hat{A}\pm 0.5\%$ of full scale
- Precise self-cleaning ball-screw
- Brush-less servo motor quarantine maintenance-free operation
- Easy calibration of load-cell
- Training video included
- Software is included
- 7/24 online distance support
- Amalgam test specimen making mold is as option (will be quote separately in case of request)



Creep Tester for Amalgam According to ISO 24234

- Amalgam test specimen making mold is as option (will be quote separately in case of request)
- Sensor for displacement measurement
- Number of stations as per customer request
- Weights for applying 36 MPa stress on sample piece
- USB port for connection to PC
- Temperature control for the chamber of sample piece
- Graphical display of temperature and displacement
- Recording data for 4 hour
- Air circulation for sample chamber
- 7/24 online distance support
- Software in included
- Windows based software
- Training video is included
- Computer is as option (Machine can work with any computer with any version of WINDOWS)

Category

1. Equipment for Standards
2. Standards